standa

S

Optical Positioners

9

Base Positioners

Translation & Rotation Stages

Fine Adjustment Screws

Opto-Mechanics & Sets

 \bigoplus

 \oplus

Extremely quiet fan <50 dBA

Odorless plastic curtains


- HEPA air filter


managed environment.


Laminar Flow Workstation 1CB is designed to provide a clean environment in any laboratory or workshop. Our filtering unit 1CB consists of a HEPA air filter and a steel frame complete with odorless plastic curtains. The vertical aligned flow of filtered air ensures that clean air flows uniformly over the workstation, keeping your samples, or optics in a clean

In a laser optics manufacturing, a laser repair or research environment, cleanliness of optical surfaces is imperative. The millions of particles present in most laboratories can fall onto optical surfaces, or they can be carried to optical surfaces by the energy in a laser beam. Once on the surface the lasing beam will burn the particles onto the coating, causing attenuation of the intra cavity energy. In addition the burnt particle will absorb more energy and cause localized heating on the surface of the optic. The localized heating can cause thermal lensing and damage to the optical coating. When adding the 1CB to your optical table you remove the particles and therefore the potential damage to your optics.

Laminar Tabletop Workstation

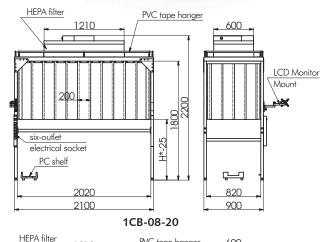

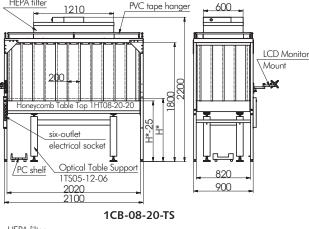

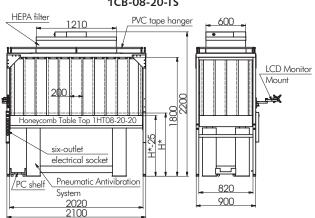

1CB-08-20-VIS

ORDERING INFORMATION

As an example model - 1CB-08-20-TS combines the proven 1HT08-20-20 (800 \times 2000 \times 200 mm) honeycomb tabletop, the sturdy 1TS05-12-06-AR table support and the NEW 1CT-08-20 air filter unit complete with curtains. For greater stability and performance we suggest ordering a Vibration Isolation System 1VIS (1CB-08-20-VIS) instead of static rigid supports 1TS.

If you already have your optical table set up and wish to add the filtered air unit to your table, we can provide any custom size to suit your needs. Quantity of HEPA filters can be increased according to the customer's application dimensions.

12

ω

9

Lasers & Laser Accessories 12

1-26

1CB-MM
Monitor mount
with articulated frame

1CB-ES
6 Europe-type
sockets

1CB-PCS
PC holder
under the table

Laminar Tabletop workstation included features:

1CB-MM – Monitor mount with articulated frame

1CB-ES - 6 Europe-type sockets1CB-WL - Day light (cool white) lamp

1CB-PCS – PC holder under the table

Featured units can be also ordered separately.

Filter specification

www.standa.lt

e-mail: sales@standa.lt

Fax: +370-5-2651483

Phone: +370-5-2651474

- Low sound, low watts, low profile, and low operating costs.
- Three speed switch features low, medium, and high settings; standard on all 2 ft \times 4 ft (600 mm \times 1210 mm) and 2 ft \times 3.5 ft (600 mm \times 1057 mm) units.
- Solid-state speed control standard on 2 ft \times 2 ft (600 mm \times 600 mm) and 2 ft \times 3 ft (600 mm \times 905 mm) units.
- Forward-inclined centrifugal-type fan.
- High Efficiency Particulate Air (HEPA) UL 900 Filter: 99.99 % efficient @ 0.3 micron.
- Extremely quiet <50 dBA.
- Snap-in pre filter allows for easy replacement and maintenance.
- Mill finished aluminum exterior.
- Tested to IEST recommended RPC standards.
- UL listed (120 V, 240 V, 277 V) with standard UL 900 filter.

ORDERING INFORMATION

Laminar Flow Workstation

Air filter unit complete with curtains, for tabletop sized 800x2000 mm.

Does not include tabletop or its legs.

Laminar Flow Workstation with rigid Table Supports

tabletop 1HT08-20-20 (800 x 2000 x 200 mm);
rigid table support 1TS05-12-06-AR;
1CB-08-20 air filter unit complete with curtains.

 $\label{lem:lemma$

pneumatic vibration isolation system 1VIS22-70-4;
 1CB-08-20 air filter unit complete with curtains.

Europe-type sockets (additional)

Monitor mount with articulated frame (additional)

PC holder under the table (additional)

Day light (cool white) lamp (additional)

1CB-08-20-TS

Includes:

1CB-08-20

1CB-08-20-VIS

Includes:

1CB-ES 1CB-MM 1CB-PCS 1CB-WL

standa



